SMARTLY[®] SAVER

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MEDICAL DEVICES

**AED
DAE
DEA**

The use of the defibrillator is indicated to treat adult and pediatric patients affected by sudden cardiac arrest. It can only be used if all of the following conditions are met:

- the patient is unconscious
- the patient is not breathing
- the patient has no heartbeat

Defibrillators can be used both in healthcare and emergency environments (ambulances, emergency rooms, etc.) and in non-healthcare environments (public or private places). The regulations governing the requirements regarding the use and training for AED devices differ from country to country. Strictly follow local laws and regulations for the use of automatic external defibrillators. In any case, it is recommended to attend a BLSD (Basic Life Support & Defibrillation) training course in order to be able to intervene effectively in an emergency both in the use of the defibrillator and in performing CPR.

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10.

Practice 2 insufflations
Voer 2 beademingen uit
Udøv 2 insufflationer

The Fully Automatic model will deliver the shock automatically, without pressing the shock button, after a warning voice message.

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
2.

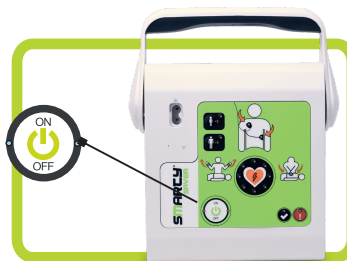
3

Start of the activation test, follow the voice messages provided by the AED

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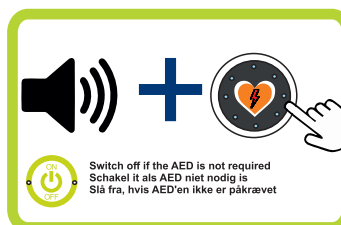Nr.	Function	Nr.	Function
1	Paediatric Mode selector To select paediatric patient to treat with universal PADS	9	Red control LED In stand-by mode: Error status
2	Adult mode selector To select adult patient to treat with universal PADS	10	PADs connector port
3	Icon "Place the PADs" Place the PADs on the patient's chest	11	SMARTY Saver microphone
4	Icon "Don't touch" Do not touch the patient	12	SMARTY Saver speaker
5	*Shock button To deliver the shock	13	Battery compartment
6	Icon CPR Start the cardiopulmonary resuscitation	14	Contact pins
7	ON/OFF button To switch ON/OFF the device	15	μSD Memory Card and USB slot
8	Green control LED In stand-by mode: Status ready use		

*In the automatic model, the shock button is replaced by the "auto shock" icon which announces the imminent automatic shock.

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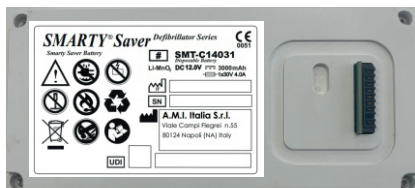
Electrodes Elektroden Elektroder

SMT-C2001

Pre-connected universal electrodes

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SMT-C2002 (FACE TO FACE)

Pre-connected universal electrodes face-to-face

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Before disconnecting the battery make sure the device is turned off for 5 seconds.

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For detailed information, always refer to the relevant user manual and to the instructions on the packaging.

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Thank you for choosing A.M.I. Italia S.r.l. **SMARTY Saver Series** defibrillator!

For the correct use of the device, it is necessary, before using it, **to read carefully the relevant user manual (downloadable from www.amiitalia.com) and this user guide.** In order to ensure a correct and rapid traceability of the product and to receive information about all updates implemented, the user is required **to register the device** in the appropriate section of our website www.amiitalia.com.

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Scan here to access the complete User Manual of A.M.I. Italia S.r.l.'s **SMARTY Saver Series** defibrillator!

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A.M.I Italia S.r.l.